
ULINE H-4737, H-4738
H-4739
OHAUS EXPLORER®
BALANCE SCALE

1-800-295-5510
uline.com

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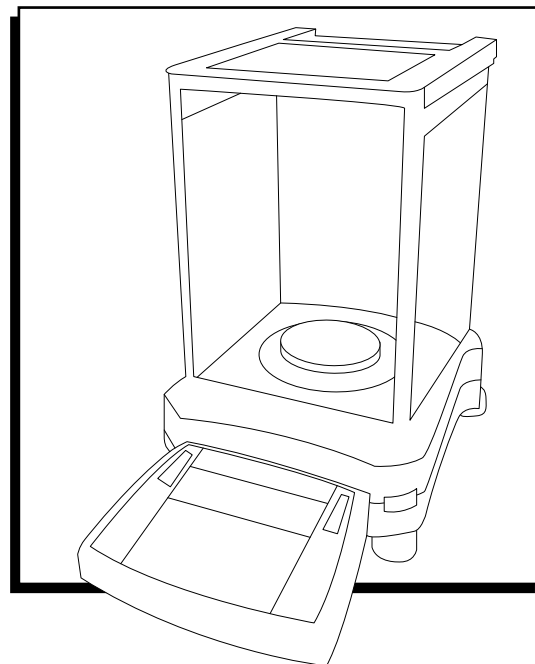OVERVIEW OF CONTROLS

CONTROLS

This balance utilizes a touch-sensitive display. Touch areas and buttons to control the functions.

OVERVIEW OF CONTROLS

FUNCTIONS AND THEIR ICONS

Seven icons are initially available to operate the balance and access the menu. For setup information, see User Settings.

Setup Icon: Touch to enter setup menu.

Standby: Touch **STANDBY** to shut down the display.

NOTE: After initial start-up, in Standby, your balance needs no warm-up time and is immediately available for weighing by touching the Start-up icon.

Printing Data: Touch **PRINT** to send the displayed value to a printer or computer via the active COM port.

NOTE: Ensure proper setup of connections, printing and interface parameters. Data may also be printed using the TouchLess sensors by configuring these for Printing. Data may also be printed using the P command from a computer connected to the COM port.

Zero Operation: Remove the load from the pan and touch **ZERO** to set the display to zero. When weighing pan is empty, the >0< indicator turns on when the measurement is within +/- 1/4 division (d) of the zero setting.

NOTE: The balance also features Auto Zero Tracking (AZT) that automatically maintains a Center of Zero within the tolerances set in the AZT menu (See Balance Settings).

TouchLess Sensor Status: Explorer balances have four TouchLess sensors that can be assigned a unique function when activated (e.g., zero, print, tare, Automatic draftshield door etc.). To assign a task for each TouchLess sensor, press **SENSORS**.

NOTE: To activate a sensor, move an object over it (from a distance of 1-2" or 3 cm). The sensor shows a green light and beeps when activated. If the sensor cannot be activated (sensor is disabled during certain instances, like when the menu is displayed), the sensor shows a red light.

AutoCal™: When AutoCal™ is on, the balance performs a self-calibration. AutoCal™ will automatically calibrate the balance (using the internal mass) each time a change in temperature significant enough to affect accuracy, or after every 11 hours. Touch **AutoCal™** to initiate. (Default status is "on".)

Menu: Touch **MENU** to enter the balance's menu list.

More...: Touch **MORE** to access Level Assist, Zero, Tare, Units, Prefare, Cal Test, Calculator and Stopwatch.

Level Assist: For assistance with leveling the balance; the Level Assist screen shows which direction to turn the thumb wheels to center the level bubble.

Applications: The balance can be configured to operate in various application modes. Touch **APPLICATIONS** to choose – or the top button in the application area.

Taring a Container: Taring sets the weight of a container as zero, so only the weight of objects held in the container (net weight) is displayed.

Manual Tare: Place the empty container on the pan and touch **TARE** when stable. Add material to the container. The net weight of the material is displayed. To clear the Tare value, remove the container from the pan and touch **TARE**.

Weighing Units: The balance can be configured to measure in a variety of weighing units, including three custom units. Touch **UNITS** to choose from the displayed list.

NOTE: The active weighing Unit can also be accessed by touching the unit area of the Main Display line; if a unit is not displayed in the list, it must first be turned on via the User Menu.

Preset Tare: To enter a preset tare, press the icon and then enter a value; enter 0 to clear the preset tare. The display will show the PT symbol and the tare value as a negative number.

Calibration Test: Used to compare a known calibration weight against the last stored span calibration data.

Calculator: To use the built-in calculator application, touch **CALCULATOR**. A four function calculator with addition, subtraction, multiplication and division mathematical operations is available.

- Perform a numeric operation using the keypad.
- Press = to view the result.
- To clear the display: Press **CE** to clear the entry, **C** to clear all.
- Touch **X** to close and return to the current application.

NOTE: Calculator is not available in software version 1.10.

OVERVIEW OF CONTROLS CONTINUED

Stopwatch: To use the built-in stopwatch application, touch **STOPWATCH**. A countdown timer with interval (lap) timer is available.

- Press **START** to begin using the stopwatch timer.
- Press **RESET** to reset the timer to zero.
- Touch **CLOSE** to close and return to the current application.

Logout: Used for user logout and return to User Login screen. This feature is only active if a user profile has been created and a user has logged in.


Save: Touch to save current setting/information.

OPERATION

MENU NAVIGATION

TO ENTER/EXIT THE MENU

All menu navigation is performed by touching the display. To enter the menu, touch **MENU** from any application home screen. The main menu appears, with buttons for **MAIN**, **BACK** and **DONE**. Continue touching the appropriate list item to navigate to the menu items or touch a button to change location.

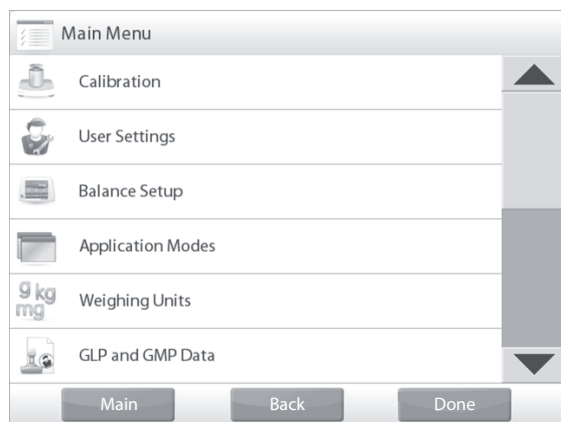CHANGE MENU SETTINGS

To change a menu setting, navigate to that setting using the following steps:

1. **Enter the Menu:** From any application screen, touch **MENU**. The main menu list appears on the display.
2. **Select Sub-Menu:** Scroll down to the selected item of the main menu list and touch it. The sub-menu appears.
3. **Select the Menu Item:** Continue until the desired setting is chosen in the menu list. Touch the setting to change it. The changed setting will be displayed as highlighted for two seconds to confirm the changed value.
4. **Exit the Menu and return to the current application:** After the setting is confirmed, touch **DONE** to return to the application.

NOTE: At any time, the **MAIN**, **BACK** and **DONE** buttons can be touched to navigate to the desired area of the menu or return to the current application.

MENU LIST

Calibration: Touch to view calibration options.

User Settings: Touch to view user preferences.

Balance Setup: Touch to view balance settings.

Application Modes: Touch to view application modes.

Weighing Units: Touch to view weighing units.

GLP and GMP Data: Insert user data for traceability.

Communication: Touch to view COM device settings and print settings.

Library: Touch to delete all library records.

Diagnostics: Touch to view diagnostic items.

Factory Reset: Touch to do a factory reset.

Lockout: Touch to access the lockout menu list.

I/O Settings: Touch to view remote I/O settings.

Information: Touch to view list of topics.

OPERATION CONTINUED

ENABLING UNITS OF MEASURE

In the balance setup menu, enter the weighing units sub menu to activate the desired units of measure.

NOTE: Due to national laws, the indicator may not include some of the units of measure listed.

mg

Milligram

g

Gram

kg

Kilogram

ct

Carat

oz

Ounce

ozt

Ounce Troy

lb

Pound

dwt

Pennyweight

Grain

Grain

N

Newton

mom

Momme

msg

Mesgal

Hkt

Tael (HK)

SGt

Tael (SG)

TWt

Tael (TW)

tical

Tical

tola

Tola

bht

Baht

C1

Custom Unit 1

C2

Custom Unit 2

C3

Custom Unit 3

TURNING A UNIT ON/OFF

1. Select the desired unit, then touch the item to toggle it either on or off.
2. The current menu item status is shown:
OFF = disabled
ON = enabled

CUSTOM UNITS

Use each custom unit to display weight in an alternative unit of measure. The custom unit is defined using a conversion factor, where the conversion factor is the number of custom units per kilogram expressed in scientific notation (Factor x 10^{Exponent}).

For example: To display weight in troy ounces (32.15075 troy ounces per kilogram), enter a factor of 0.321508 and an exponent of 2.

Set the status:

OFF = disabled

ON = enabled

When custom unit is set ON, the factor, exponent and least significant digit must be set.

FACTOR

Set the conversion factor using the numeric keypad:

Settings of 0.00001 to 1.9999999 are available. The default setting is 1.000000.

EXPONENT

Set the factor multiplier:

0 = multiply the factor by 1 (1 x 10⁰)

1 = multiply the factor by 10 (1 x 10¹)

2 = multiply the factor by 100 (1 x 10²)

3 = multiply the factor by 1,000 (1 x 10³)

-1 = divide the factor by 10 (1 x 10⁻¹)

-2 = divide the factor by 100 (1 x 10⁻²)

-3 = divide the factor by 1,000 (1 x 10⁻³)

LEAST SIGNIFICANT DIGIT

Set the graduation:

Settings of 0.00001, 0.00002, 0.00005, 0.0001, 0.0002, 0.0005, 0.001, 0.002, 0.005, 0.01, 0.02, 0.05, 0.1, 0.2, 0.5, 1, 2, 5, 10, 20, 50, 100, 200, 500 and 1,000 are available.

NOTE: Least significant digit setting selections are dependent on the factor and exponent settings. Selections are limited.

CALIBRATION

The balance offers a choice of six calibration methods: Internal Calibration, Automatic Calibration (AutoCal), Calibration Adjust, Span Calibration, User Calibration and Calibration Test.

 **NOTE:** Do not disturb the balance during calibration.

INITIAL CALIBRATION

When the balance is first installed, and when it is moved to another location, it must be calibrated to ensure accurate weighing results. The balance has built-in autocal, which can calibrate the balance automatically and does not require calibration masses. If preferred, the balance can be manually calibrated with external masses. Have the appropriate calibration masses available before beginning calibration.

AUTOCAL™ INTERNAL CALIBRATION

Calibration is accomplished with the internal calibration mass. Internal calibration can be performed at any time, provided the balance has warmed up to operating temperature and is level.

With the balance turned ON and no load on the pan, touch **AUTOCAL™ INTERNAL CALIBRATION** to begin calibration.

Alternately, touch the shortcut **CALIBRATION** icon from any Application home screen to initiate an Internal Calibration.

The display shows the status, then returns to the current application. To cancel at any time, press **CANCEL**.

AUTOMATIC CALIBRATION

When Automatic Calibration is on, the balance performs a self calibration when it senses a temperature change of 34.7° F or every 11 hours.

AutoCal will automatically calibrate the balance (using the internal mass) each time there is a change in temperature significant enough to affect accuracy.

AUTOCAL ADJUSTMENT

Use this calibration method to adjust the zero calibration point, without affecting the span or linearity calibration.

Calibration Adjust may be used to adjust the result of the internal calibration by +/- 100 divisions.

NOTE: Before making a calibration adjustment, perform an internal calibration. To verify whether an adjustment is needed, place a test mass equal to the span calibration value on the pan and note the difference (in divisions) between the nominal mass value and the actual balance reading. If the difference is within +/- 1 division, calibration adjustment is not required. If the difference exceeds +/- 1 division, calibration adjustment is recommended.

Example:

Actual Weight Reading	200.014
Expected Weight Reading	200.00 (Test mass value)
Difference Weight (d)	0.014
Difference Weight in Digits	-14 (Adjust value)

To perform a Calibration Adjustment, touch **AUTOCAL ADJUSTMENT** from the Calibration Menu. Enter the value (positive or negative divisions) to match the difference noted earlier in the procedure and press **-**.

Recalibrate using internal calibration. After calibration, place the test mass on the pan and verify that the reading agrees with the test mass. If not, repeat steps until the reading matches the test mass.

Once completed, the balance stores the adjustment value and the display returns to the current application.

SPAN CALIBRATION

Span calibration uses two calibration points, one at zero load and the other at specified full load (span).

CALIBRATION WEIGHTS		
ULINE MODEL NO.	OHAUS MODEL NO.	SPAN CAL POINTS (G)
H-4737	EX224	100g, 150g, 300g
H-4738	EX4202	2,000g, 3,000g, 4,000g
H-4739	EX6202	3,000g, 4,000g, 5,000g , 6,000g

NOTE: Factory default settings are shown in bold.

With the balance turned on and no load on the pan, touch **SPAN CALIBRATION** to initiate the procedure. Additional calibration values to be used are shown on the display. The best accuracy is achieved using the mass closest to the full span value.

To change the span calibration point, touch the alternate weight shown on the display. Follow the screen instructions and place the specified calibration weight on the balance when prompted to do so. When complete, the display shows the span calibration status and returns to the current application.

USER CALIBRATION

Use this calibration method to calibrate the balance using a defined mass.

1. With no load on the scale, press **USER CALIBRATION** from the calibration sub-menu. The display shows the current full load (span) calibration point and calibration unit of measure.
2. To change the value to a user-defined weight, touch **ENTER VALUE**. Enter the weight via the pop-up screen. The balance then prompts for the user-defined weight.
3. Place the specified calibration weight on the balance and the calibration process completes. The display shows status, then returns to the current application.


CALIBRATION TEST

Use this calibration method to compare a known calibration weight against the stored span calibration data.

1. With no load on the balance, press **Cal Test** to begin the process.
2. The balance captures the zero point, then prompts for the span value.
3. The display shows status, followed by the difference between the current calibration weight and the stored calibration data.
4. After three seconds, the display returns to the current application.

TURNING AN APPLICATION ON/OFF


1. In the balance Setup Menu, enter the Application Modes sub-menu to activate the desired balance applications. Only one can be running at a time.
2. Touch the line of the desired application.

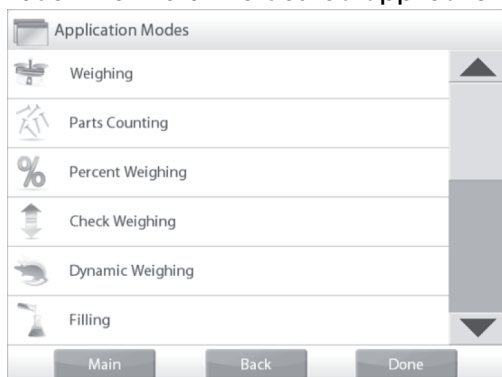3. In the item option screen, touch **ENABLED** to toggle it on or off.

4. Once an application is enabled (turned on), it may be chosen from the Applications shortcut (top-left of all home screens).
5. The current menu item status is shown:
OFF = disabled
ON = enabled


APPLICATIONS

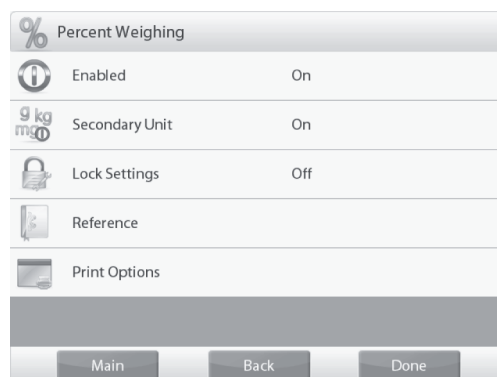The balance can be configured to operate in various application modes. Touch **APPLICATIONS** to choose – or the top button in the application area.

NOTE: Before using any application, ensure the balance has been leveled and calibrated. The balance is shipped with all units enabled.

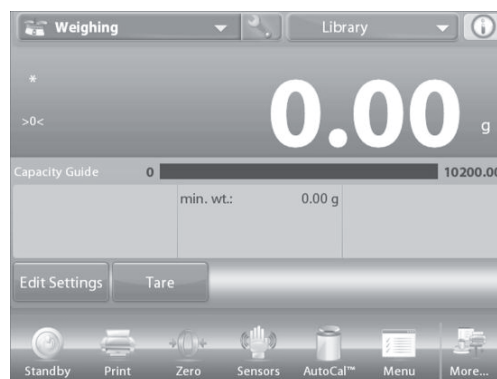WEIGHING

Use this application to determine the weight of items in the selected unit of measure.


1. In the upper left portion of the home screen, select **WEIGHING** (this application is the default).
2. Press **TARE** or **ZERO** if necessary to begin.
3. Place objects on the pan to display the weight. When stable, the * appears.
4. The resulting value is displayed in the main weighing line in the active unit of measure.

EDIT SETTINGS

To view or adjust the current settings:

1. Touch the **EDIT SETTINGS** button. The settings screen appears.
2. Minimum Weight: Establish a minimum weight value to be used to verify a reading. If an actual weight is below the established Minimum Weight value, it is flagged by a color change: yellow.
3. To adjust the Minimum Weight value, touch the **MINIMUM WEIGHT** button.

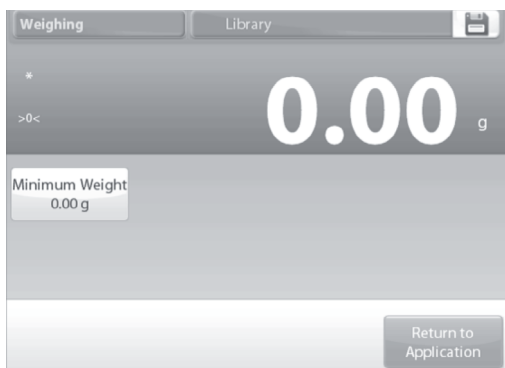4. A numeric input window appears. Use the keys to enter the desired Minimum Weight, then press $\sqrt{}$.
5. The display reverts to the previous screen with the new Minimum Weight value highlighted in orange.
6. To return to the Weighing home screen, touch **RETURN TO APPLICATION** at the bottom-right of the screen.

PARTS COUNTING

Use this application to count samples of uniform weight. Three different modes can be selected: Standard, Check, or Fill Counting.

STANDARD PARTS COUNTING (DEFAULT)

1. In the upper left portion of the home screen, select **PARTS COUNTING**. The default (or last) average piece weight (APW) is displayed.
2. Place objects on the pan to display the number of pieces.


NOTE: Auto Optimization improves counting accuracy. It recalculates the piece weight automatically as parts are added. By default, auto optimization is set to "off".

EDIT SETTINGS

To view or adjust the current settings, touch the **EDIT SETTINGS** button. The settings screen appears:

- Settings available: **APW, SAMPLE SIZE**
- Functions available: **RECALCULATE APW, RETURN TO APPLICATION**

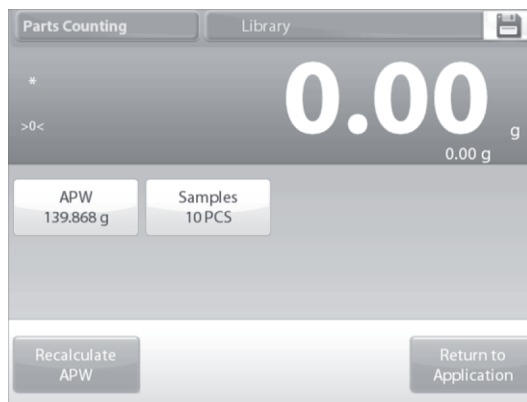OPERATION CONTINUED

ESTABLISH AN AVERAGE PIECE WEIGHT (APW)

Each time a new type of part is counted, the nominal weight of one piece (average piece weight or APW) must be established using a small quantity of pieces. This APW is stored until replaced by another APW.

Two methods establish the APW value:

- The actual APW is known.
- The APW must be derived by weight. For this case the current sample size will be used.

To choose a different value, establish it beforehand, using the **SAMPLE SIZE** button. (See Set a Known Average Piece Weight.)

SET A KNOWN AVERAGE PIECE WEIGHT (APW)

1. To adjust the APW value directly, touch the **APW** button. A numeric input window appears.
2. Key in the desired APW, then press $\checkmark$. The display returns to the previous screen with the new APW value highlighted.


3. To return to the counting home screen, touch **RETURN TO APPLICATION**.


SET A NEW AVERAGE PIECE WEIGHT (APW) – DERIVED

1. To establish a new APW, touch the **RECALCULATE APW** button. The APW Screen is displayed.

 **NOTE:** The sample size that is displayed will be used. To use a different sample size, change that first.

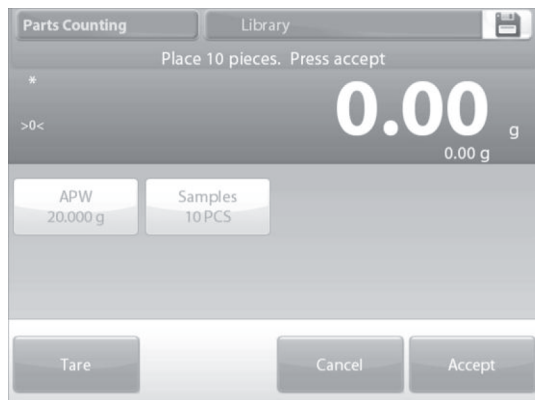2. Follow the screen instructions to establish a new APW. Actions available: **TARE, CANCEL OR ACCEPT**.
3. Press **ACCEPT** to capture the value and return to the settings screen. The settings screen returns with the new APW value highlighted.

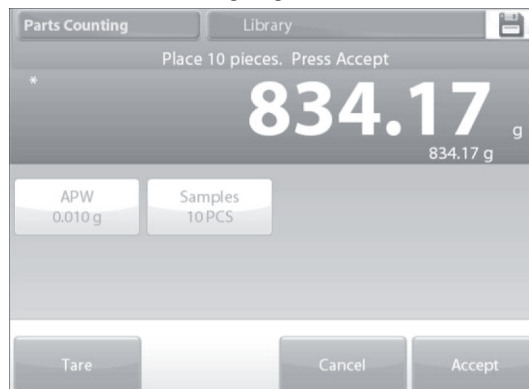4. To return to the counting home screen, touch **RETURN TO APPLICATION**. The home screen shows 10 pieces at the new APW.


SAMPLE SIZE

The sample size can be 1 to 10,000 pieces. The default sample size is 10. Once a sample size is changed, the balance will immediately open the recalculate APW screen, expecting to establish a new APW.


The current sample size is shown on the settings screen.

1. To adjust the sample size, touch the **SAMPLES** button. A numeric input window appears.

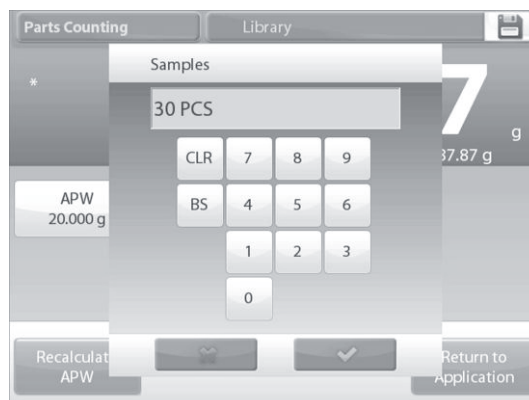OPERATION CONTINUED

2. Key in the desired sample size, then press $\checkmark$.
3. The next screen appears, with the message to place the number of pieces entered in the numeric input window.
4. Place the indicated number of pieces on the platform. Touch either **TARE**, **CANCEL** or **ACCEPT**.

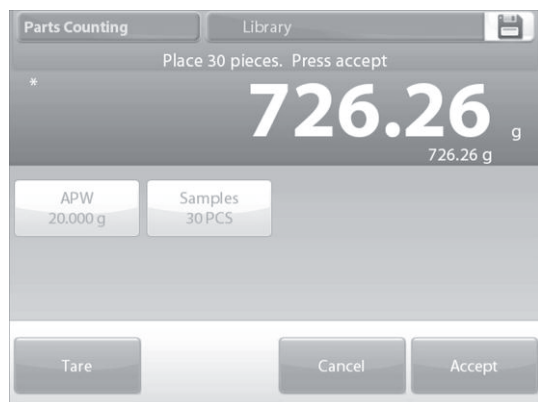
5. Press **ACCEPT** to capture the value and return to the settings screen.
6. The settings screen appears with the new APW value highlighted.
7. To begin counting pieces, touch **RETURN TO APPLICATION**.


PARTS COUNTING – CHECK

This feature permits establishing a set quantity of pieces as criteria for similar items that can be quickly checked against a sample.


1. Press the **SETUP** button to enter the counting setup menu.

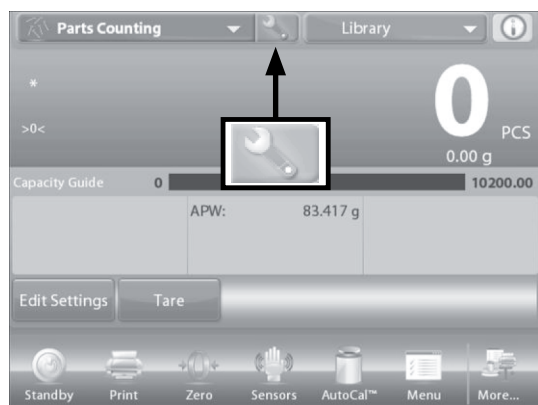2. Change counting mode to check counting.

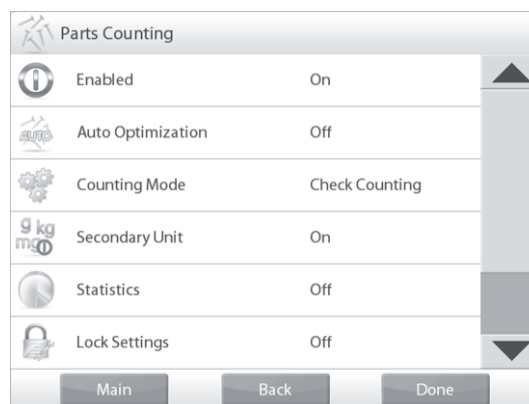
3. The default (or last) APW is displayed. Place objects on the pan to display the check status (number of pieces).

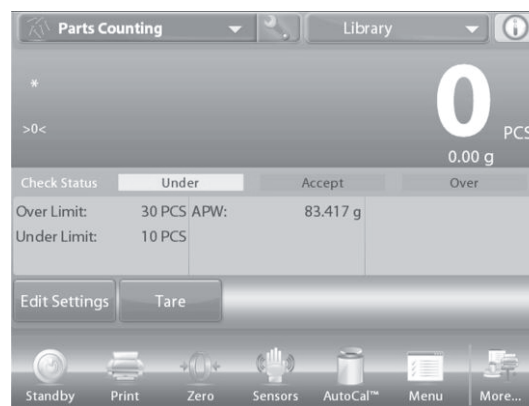EDIT SETTINGS

To use the parts counting-check mode, the over limit and under limit values must be set. The APW and standard sample size can also be set.

1. Touch the **EDIT SETTINGS** button. The settings screen appears. The buttons for **APW**, **SAMPLES**, **OVER LIMIT** and **UNDER LIMIT** show default (or previously set) values. To change any value, touch the button. A keyboard screen appears.

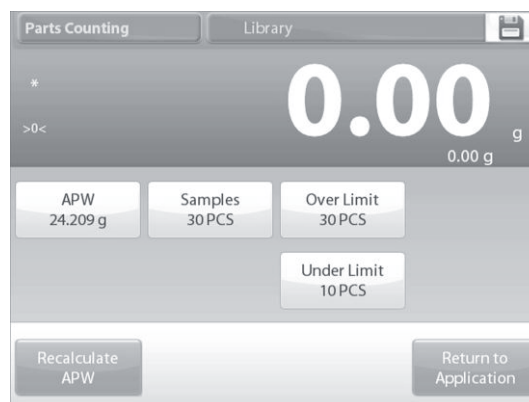2. To set a new APW, press the **RECALCULATE APW** button.

OPERATION CONTINUED

3. A new screen appears with the message: Place 30 pieces. Press **ACCEPT**.

NOTE: Default is 10 pieces, but message matches number of pieces in **SAMPLES** button.


4. After following these instructions, press **RETURN TO APPLICATION** on the main edit settings screen.

PARTS COUNTING – FILL

This feature permits establishing a set quantity of pieces as a reference, so when pieces are added, their percentage of the reference amount appears. A Progress bar indicates the percentage.


1. In the upper left portion of the home screen, select **PARTS COUNTING**. (The option will appear if it has been enabled in the application modes menu.)
2. Press the **SETUP** button to enter the counting setup menu.
3. Change counting mode to fill counting.
4. The default (or last) APW is displayed. Place objects on the pan to display the fill status.

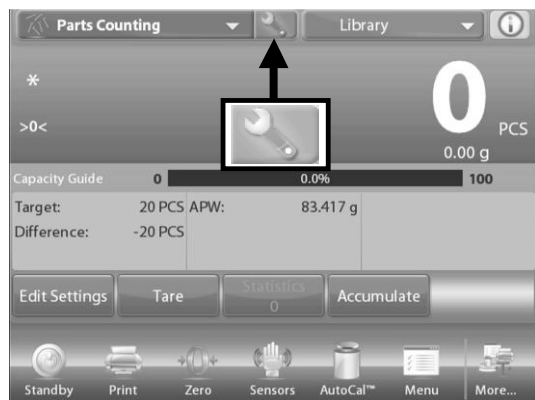
EDIT SETTINGS

Touch the **EDIT SETTINGS** button. The settings screen appears.

The buttons for **APW**, **SAMPLES** and **TARGET** show default (or previously set) values.

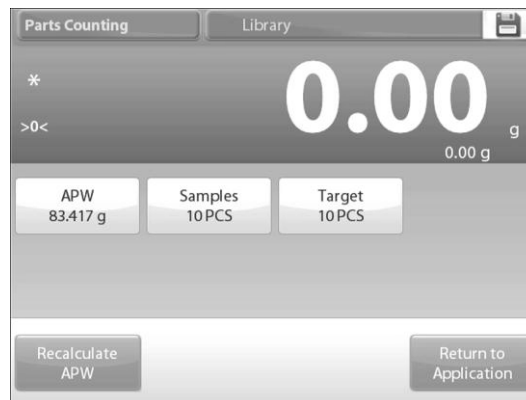1. To change values, touch the button. A keyboard screen appears.
2. Key in the new value, and press $\checkmark$.
3. To set a new average piece weight, press **RECALCULATE APW**.
4. A new screen appears saying: Place reference weight on the pan. Press **ACCEPT** to continue.
5. This establishes a new reference weight with a sample. Alternatively, fill in values in the screens for APW, samples and target.

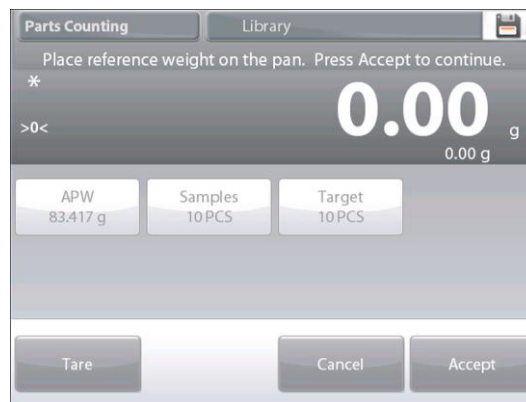6. After following these instructions, and setting all values as desired, press **RETURN TO APPLICATION**.

PERCENT WEIGHING

Use percent weighing to measure the weight of a sample displayed as a percentage of a pre-established reference weight.


1. The default (or last) reference weight is displayed. Place an object on the pan. The difference between the sample and the reference weight is displayed as both a weight and a percentage.
2. To access percent weighing press the upper left portion of the home screen, so the available application modes are visible.
3. Select **PERCENT WEIGHING**. The percent weighing home screen appears.

EDIT SETTINGS

To view or adjust the current settings, touch the **EDIT SETTINGS** button. The settings screen appears:

- Settings available: **REFERENCE WEIGHT, REFERENCE ADJUST**
- Functions available: **RECALCULATE REFERENCE, RETURN TO APPLICATION**

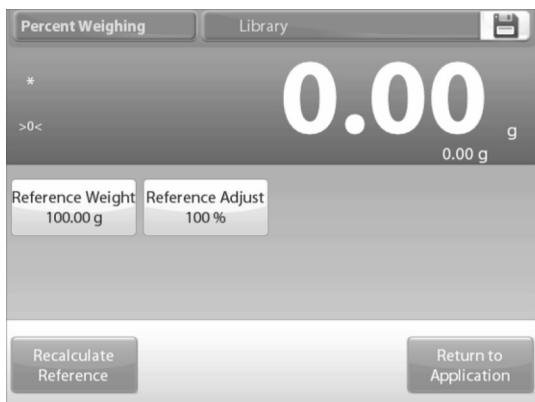ESTABLISHING A REFERENCE WEIGHT

There are two methods to establish the Reference Weight value:

- **Direct:** If the actual Reference Weight is known, it can be entered directly.
- **Derived:** The sample's weight can be used to establish the reference weight.

ESTABLISHING A KNOWN REFERENCE WEIGHT

1. The current Reference Weight value is shown on the Settings screen. To adjust the Reference Weight directly, touch the **REFERENCE** button. A numeric input window appears.
2. Key in the desired Reference Weight, then press $\sqrt{}$. The display returns to the previous screen with the new value highlighted.
3. To return to the Percent Weighing home screen, touch **RETURN TO APPLICATION**.

USING A SAMPLE TO ESTABLISH A DERIVED REFERENCE WEIGHT

1. To establish a new Reference Weight value, touch the **RECALCULATE REFERENCE** button.
2. The New Reference screen appears. Follow the screen instructions to establish a new Reference Weight. Actions available: **TARE, CANCEL** or **ACCEPT**.

3. Press **ACCEPT** to capture the value and return to the Settings screen.
4. The display returns to the Settings screen with the new Reference weight value highlighted momentarily.
5. To return to the Percent weighing home screen, touch **RETURN TO APPLICATION**.

OPERATION CONTINUED

REFERENCE WEIGHT ADJUSTMENT (REFERENCE ADJUST)

The Reference Weight adjustment can be 1 to 1,000%.

Examples: 100% means the sample is displayed as compared to 100% of the Reference Weight.
25% means the sample is displayed as compared to 25% of the Reference Weight.

The Edit Settings screen shows current Reference Adjust value. To change it, touch **REFERENCE ADJUST** button.

1. A numeric window appears. Key in the desired Reference Weight Adjustment, then press $\sqrt{}$.
2. The display returns to the previous screen with the new value highlighted momentarily. To return to the Percent Weighing home screen, touch **RETURN TO APPLICATION**.

CHECK WEIGHING

Check weighing is used to compare the weight of a sample against target limits. Three different modes can be selected: Over and Under, Nominal Weight Tolerance or Nominal Percent Tolerance.

STANDARD CHECK WEIGHING (DEFAULT)


1. In the upper left portion of the home screen, select **CHECK WEIGHING**.
2. The default (or last) check weigh limits are displayed. Place objects on the pan. The Under/Accept/Over status is shown in the progress bar area, while the actual weight of the item is shown on the main display line.

EDIT SETTINGS

To view or adjust the current settings touch the **EDIT SETTINGS** button. The Settings screen appears:

- Settings available: **OVER LIMIT, UNDER LIMIT**
- Functions available: **RETURN TO APPLICATION**

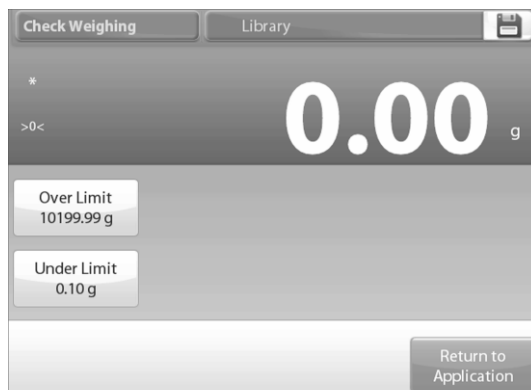3. The check limit values are displayed in the settings screen.
 - To set the over limit value, touch the **OVER LIMIT** button.
 - To set the under limit value, touch the **UNDER LIMIT** button.
4. A numeric input window appears. Enter the desired limit weight, then press $\sqrt{}$.


5. The display returns to the previous screen with the new value highlighted momentarily. To return to the check weighing home screen, touch **RETURN TO APPLICATION**.

OPERATION CONTINUED

CHECK WEIGHING – NOMINAL/WEIGHT TOLERANCE MODE

With this feature, you can establish a Target or Nominal weight value and pre-set tolerance in weight.

1. To enable this feature, enter the setup menu.
2. Change the mode to Nominal, wt. Tolerance. The Check Weighing options screen returns.

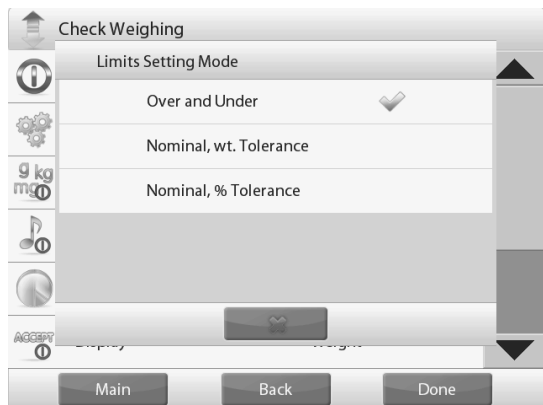3. Press **DONE** to begin check weighing.
4. To use check weighing, touch **EDIT SETTINGS** to set the nominal, and upper and lower tolerance values (in weight).
5. Press **RETURN TO APPLICATION** when settings are complete.

CHECK WEIGHING – NOMINAL/PERCENT TOLERANCE MODE

With this feature you can establish a target or nominal weight value and pre-set tolerance in percent.

1. To enable this feature, enter the setup menu.
2. Change the mode to Nominal, % Tolerance. The Check Weighing options screen returns. Press **DONE** to begin check weighing.

3. To use this mode, touch **EDIT SETTINGS** to set the target weight, and upper and lower tolerance values.
4. Press **RETURN TO APPLICATION** when settings are complete.

DYNAMIC/ANIMAL WEIGHING

Use this application to weigh an unstable load, such as a moving animal. Three different start/reset modes can be selected: manual (start and stop via key press), semi-automatic (auto-start with manual reset), and automatic (start and stop automatically).


EDIT SETTINGS

To view or adjust the current settings, touch the **EDIT SETTINGS** button. The Application Settings screen appears:

- Settings available: **AVERAGING TIME**
- Functions available: **RETURN TO APPLICATION**

1. The default value is 10. To change it, press the button.


 **NOTE:** The Averaging Time can be set from 1 to 99 seconds.

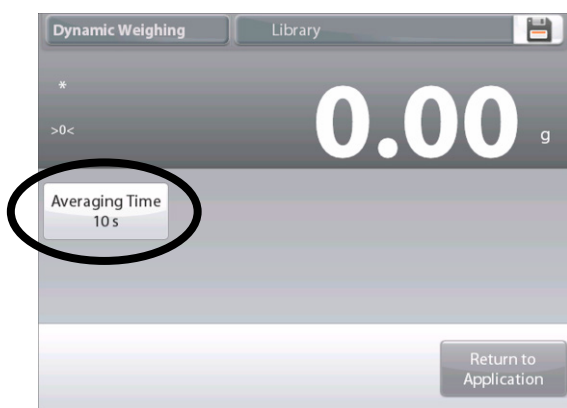
2. A numeric input window is displayed. Key in the desired Averaging Time, then press **DONE**.
3. The display returns to the previous screen with the new value highlighted. To return to the Dynamic Weighing home screen, touch **RETURN TO APPLICATION**.

DYNAMIC WEIGHING – MANUAL (DEFAULT)

1. To set the mode, touch the **SETUP** wrench. When you make a selection, the setup options screen reappears. Press **DONE**.

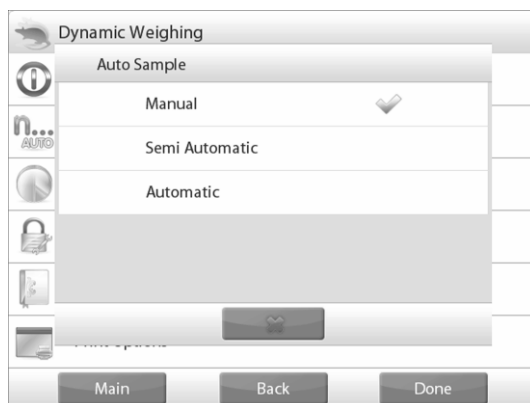2. On the Dynamic Weighing home screen, the default (or last) averaging time is displayed. Follow the screen instructions.

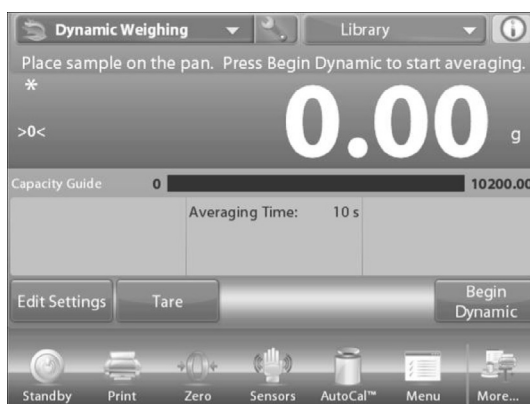3. Place objects on the pan. Press the **BEGIN DYNAMIC** button.
4. The balance begins a countdown (averaging process). During the countdown, the information line displays the time remaining. If necessary, press **CANCEL** to quit and return to the home screen.

OPERATION CONTINUED

- When the countdown ends, the result is displayed and held.

- Press **RESET** to clear the held weight and return to the initial screen.

NOTE: Clear the pan before beginning a new Dynamic Weighing cycle.

DYNAMIC WEIGHING – SEMI-AUTO

The cycle begins when an object (animal) is placed on the pan, but the held value must be reset manually.

- To enable this feature, touch the **SETUP** icon in the top-center of the main screen to enter the Dynamic Weighing setup menu.


- Change the mode to Semi Automatic.

- The dynamic – semi-automatic screen appears with the default (or last) averaging time displayed.

- Place specimen on the pan. The cycle starts automatically.

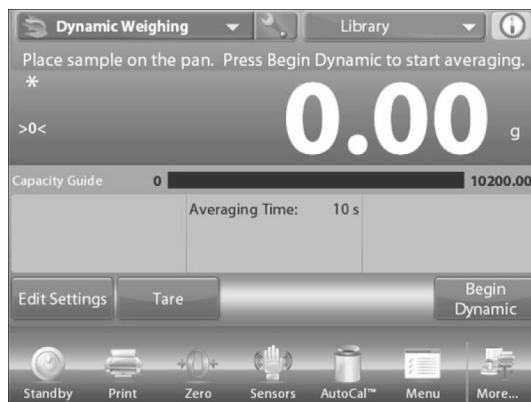
- After the cycle ends, the weight reading is held. Press the **RESET** button to clear the reading and return to zero.

DYNAMIC WEIGHING – AUTOMATIC


The cycle begins when an object (animal) is placed on the pan, and the held value is automatically reset when the object is removed from the pan (10-second delay).

- To enable this feature, enter the Dynamic Weighing setup menu.
- Press **AUTO SAMPLE**, and change the mode to Automatic.

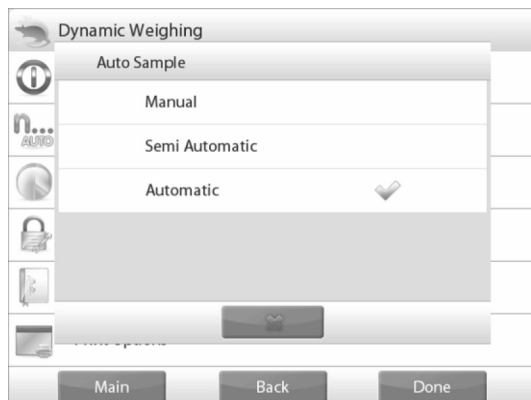OPERATION CONTINUED

- The dynamic – automatic screen appears with the default (or last) averaging time displayed. Place objects on the pan to begin. The cycle starts automatically.

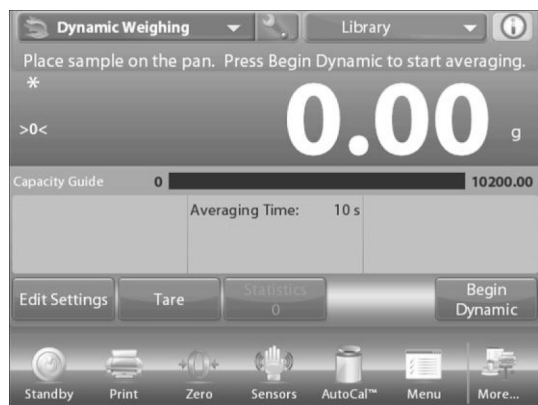- After the cycle ends, remove the specimen. The balance returns to the home screen automatically (after a 10-second delay).

- Press the **RESET** button to clear the reading and return to zero.

FILLING

This application allows the user to slowly fill a container to a pre-determined target weight. The progress bar displays the filling status, and within 10 percent of the target value the progress bar converts to fine resolution (+/- 10%) for accurate results.

- In the upper left portion of the home screen, select **FILLING**.

- The default (or last) target weight is displayed. Place objects on the pan to begin. The progress bar displays the filling status, and within 10 percent of the target value, the progress bar converts to fine resolution display.

EDIT SETTINGS

To view or adjust the current settings, touch the **EDIT SETTINGS** button. The settings screen appears:

- Settings available: **TARGET WEIGHT**
- Functions Available: **RECALCULATE TARGET WEIGHT, RETURN TO APPLICATION**

ESTABLISH A TARGET WEIGHT

Each time a new item is used for filling, the Target Weight (fill-to-weight) must be established.


There are two methods to establish the Target Weight:

- The actual Target Weight is known.
- The Target Weight can be derived from a sample.

OPERATION CONTINUED

ESTABLISH USING A KNOWN TARGET WEIGHT

The current Target Weight appears on the **Settings** screen, on the **TARGET WEIGHT** button.

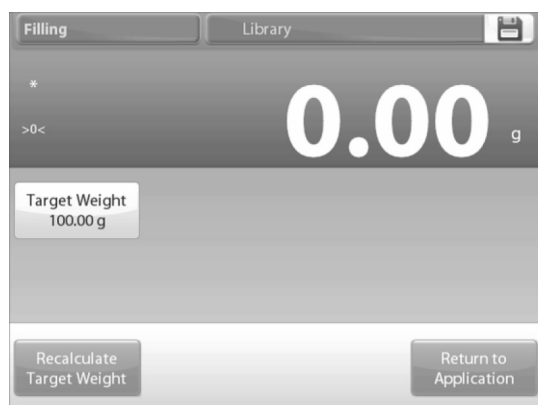1. To adjust the Target Weight value directly, touch the **TARGET WEIGHT** button.
2. A numeric input window appears. Enter the desired Target Weight, then press $\checkmark$. The display returns to the previous screen with the new value highlighted.
3. Press **RETURN TO APPLICATION** to go to the Filling home screen.

USING A SAMPLE TO ESTABLISH A NEW TARGET WEIGHT – DERIVED

1. To establish a new Target Weight value, press the **RECALCULATE TARGET WEIGHT** button.
2. The New Target Screen appears, with a message: "Place target weight on the pan". Press **ACCEPT** to continue.

3. Place the weight on the pan (or in a tared container).
4. Press **ACCEPT** to return to the Settings screen.
5. The display returns to the Settings screen with the new Target Weight value highlighted momentarily.
6. To return to the Filling home screen, touch **RETURN TO APPLICATION**.

TOTALIZATION

Totalization measures the cumulative weight of a sequence of items. Two different totalization accumulation modes can be selected: manual accumulation (requires key press) and automatic accumulation. The cumulative total may exceed the capacity of the balance.

TOTALIZATION – MANUAL (DEFAULT)

1. Touch the upper left portion of the home screen.
2. A drop-down menu appears. Select **TOTALIZATION** if it is present. (If not, touch the Menu icon in the bottom-left corner, then choose **APPLICATION MODES**, scroll to **TOTALIZATION**, press it, then press **ENABLED** to turn it on. It will now appear in the application modes menu on the home screen).

3. Place objects on the pan to begin.
4. The sample weight is shown on the main display line. Press **ACCUMULATE** to add the weight of the item to the total.


5. When stable, the new total is shown on the secondary display line.
6. Remove the item from the weighing pan, then add the next item and continue as above.
7. When finished, press the **CLEAR TOTAL** button to reset the second display line to zero.

OPERATION CONTINUED

TOTALIZATION – AUTOMATIC

The weight of the item is added automatically when stable.

1. To enable this feature: Touch the **SETUP** icon to enter application setup.
2. Press **AUTO SAMPLE** to turn this feature on.
3. Clear the pan if necessary. Add a new sample. The item weight is shown on the main display line.

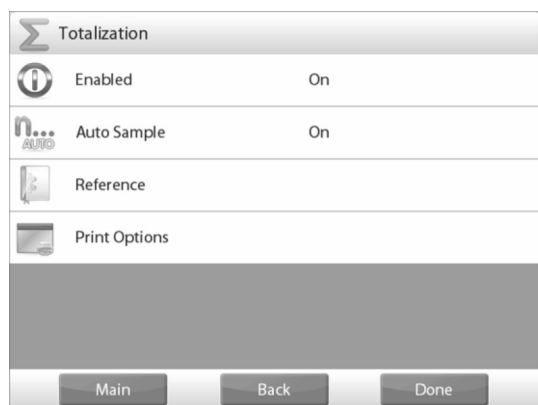4. The weight of the item is automatically added, when stable, to the total. The new total is shown on the second display line.
5. Remove the item from the weighing pan, then add the next item and continue the same as in manual mode.
6. When finished, press the **CLEAR TOTAL** button to reset the second display line to zero.

FORMULATION

Use this application for compounding and recipe making. The number of components can be 2 to 99. Formulation has two available modes of operation: free formulation and recipe formulation.

FREE FORMULATION (DEFAULT)

This mode allows the user to freely add components. A filler material is asked for at the end.

1. In the upper left portion of the home screen, select **FORMULATION**.
2. To start the compounding process, press **BEGIN FORMULATION**.


3. A prompt says place component 1 then press **ACCEPT**. The item weight is displayed on the main display line. Place the first ingredient on the pan (or in a tared container). Press **ACCEPT** to store the first component.

4. A prompt says place item 2. Press **ACCEPT**. Continue adding components and pressing **ACCEPT** to store the weight of the individual components until the formula is complete.
5. The total line shows the total weight of all the components. Press **PRINT** at any time in the process to print the recipe.
6. Then press **END FORMULA**. The display may prompt for a filler material. If the formula calls for a filler material, press **OK**.

OPERATION CONTINUED

7. A pop-up message window appears, prompting for the filler material to be added. Add the filler if needed, then press **ACCEPT**. (If no filler is needed, press **ACCEPT** without adding filler.) The results are then displayed in a pop-up window.

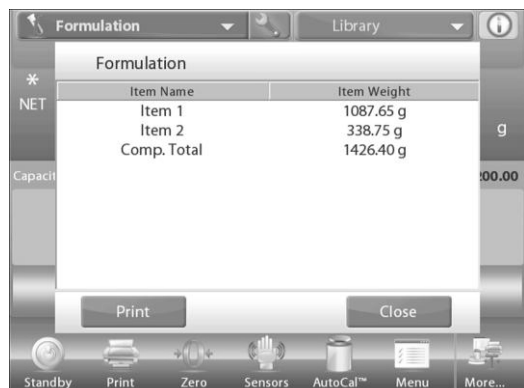
FREE FORMULATION RESULTS

1. Press **PRINT** to send data to the enabled COM port.
2. Press **CLOSE** to close the results window and return to the application home screen.

RECIPE FORMULATION


Recipe formulation requires the formula to be entered first, followed component by component.

1. To enable this feature, press the **SETUP** icon to enter the formulation setup menu.
2. Change the mode to Recipe Based.

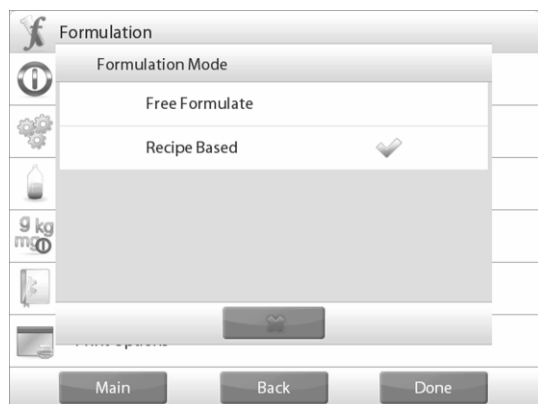
3. Press **DONE** in the main setup options screen.
4. Press **BEGIN FORMULATION**.

5. If no items have been previously set up, the settings screen appears. If items already exist in memory, formulation begins. Press the **ADD NEW** button.

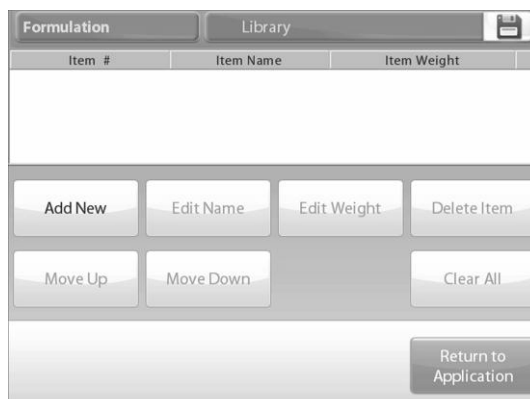6. A new default component is added to the end of the formula (default name, weight). To change the default component name, press **EDIT NAME**. The alpha-numeric entry screen appears.

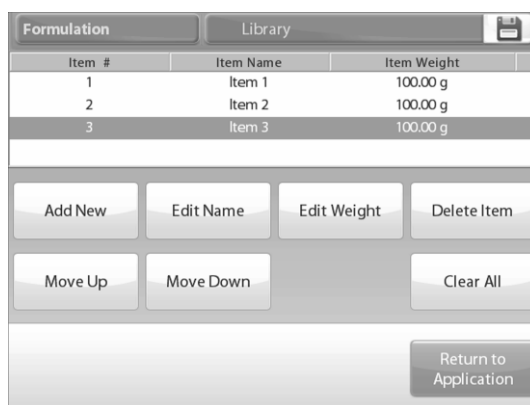7. Enter the new name, then press **√**. The display returns to the previous screen with the new value highlighted momentarily.
8. To change the default component weight, press **EDIT WEIGHT**. The numeric entry screen appears.
9. Enter the new weight, then press **√**.
10. The display returns to the previous screen with the new value highlighted momentarily.

OPERATION CONTINUED

- Continue adding components as required for your formula.
- To change the order of a component in the formula, touch the item to highlight it, then press **MOVE UP** or **MOVE DOWN**.

Item #	Item Name	Item Weight
1	Item 1	100.00 g
2	Item 2	100.00 g
3	calcium	100.00 g

Buttons: Add New, Edit Name, Edit Weight, Delete Item, Move Up, Move Down, Clear All, Return to Application

- To delete a component, touch the item to highlight it, and press **DELETE ITEM**.
- To clear all components, press **CLEAR ALL**. A confirmation is required.
- When finished, press **RETURN TO APPLICATION**.
- Press **BEGIN FORMULATION**. Follow the screen instructions: Add the required component in the formula.

Step 1: Place 20.00 g Item 1 Press Accept.

499.33 g

Fill Status: 90 (2496.7%) 110

Diff. (wt.): 479.33 g Target: 20.00 g

Diff. (%): 2396.7 %

Buttons: Item Settings, Tare, Cancel, Accept

Bottom Bar: Standby, Print, Zero, Sensors, AutoCal™, Menu, More...

- The component weight is shown on the main display line. The progress bar displays the fill status. Within 10% of the target value the progress bar converts to fine resolution display.
- Press **ACCEPT** to add the weight of the component to the total. The new total appears on the second display line.
- The difference to the target is displayed for each component (in weight and percent). Follow the screen instructions to continue adding components and pressing **ACCEPT** to complete the formula.
- When finished, the formulation results screen is displayed.

FORMULATION RESULTS

- Press **PRINT** to send the results to the enabled COM port.
- To close the results window and return to the Application home screen, press **CLOSE**.

Item Name	Target	Actual	Diff. (%)
Item 1	20.00 g	499.33 g	2396.7 %
Item 2	100.00 g	784.61 g	684.6 %
calcium	100.00 g	820.13 g	720.1 %
Comp. Total	220.00 g	2104.07 g	856.4 %

Buttons: Print, Close

Bottom Bar: Standby, Print, Zero, Sensors, AutoCal™, Menu, More...

DIFFERENTIAL WEIGHING

Differential weighing stores weight values so samples can be dried or processed and the difference in weight calculated. Up to 99 samples can be stored.

Differential weighing has two operating modes:

- Auto-Process ON, where the initial and final weights are sequentially gathered for a set of items, with no need to return to item settings.
- Auto-Process OFF, where initial and final weights are gathered separately for each item, after each item is selected in item settings.

Differential

Enabled On

Auto-Process Off

Reference

Print Options

Buttons: Main, Back, Done

Touch the **SETUP** icon to enter application setup from the home screen.

ITEM SETTINGS

- Press **ITEM SETTINGS** to set up items in memory.
- To set up the item, press **ADD NEW**.

OPERATION CONTINUED

3. New items are added with default name and weight information (no initial or final weights).
4. To change the default item name, press **EDIT NAME**. The alpha-numeric entry screen appears.

Item Name	Initial wt.	Final wt.	Difference
Item 1	---	---	---
Item 2	---	---	---
Item 3	---	---	---

Buttons: Edit Name, Clear Initial Weight, Clear Final Weight, Delete Item, Add New, Clear All, Return to Application

5. Enter the new name, then press $\checkmark$. The display returns to the previous screen with the new value highlighted momentarily.
6. Continue adding items as required using **ADD NEW**. New items can be added at any time.
7. To delete an item, touch the item to highlight it, then press **DELETE ITEM**.
8. To clear the item initial weight press **CLEAR INITIAL WEIGHT**. This feature is disabled if no initial weight is stored.
9. To clear the item's final weight, press **CLEAR FINAL WEIGHT**. This feature is disabled if no final weight is stored.
10. To clear all items, press **CLEAR ALL**, then confirm.
11. When finished, press **RETURN TO APPLICATION**.
12. On the main Application screen, follow the screen instructions. Add the first item, then press **ACCEPT** to store the value.
13. The initial weights of the defined item are recorded first, then the final weight so the items are collected next.
14. Once the item's initial weight and the final weight are stored, the values are displayed along with the difference in weight and percent.
15. The items are stored in memory. To begin a new Differential Weighing session, delete the initial and/or final weights.

DIFFERENTIAL WEIGHING WITH AUTO-PROCESS ON

Differential weighing with auto process ON requires the initial and final weights to be sequentially gathered for a single item.

1. To enable this feature, touch the **SETUP** icon from the home screen to enter application setup.
2. Press **AUTO PROCESS** to turn it on.
3. Follow the screen instructions to begin.
4. Press **ITEM SETTINGS** to set up items in memory if necessary. Use the buttons as necessary.

Item Name	Initial wt.	Final wt.	Difference
Item 1	---	---	---
Item 2	---	---	---
Item 3	---	---	---

Buttons: Edit Name, Clear Initial Weight, Clear Final Weight, Delete Item, Add New, Clear All, Return to Application

5. Follow the steps as shown on page 20. Press **RETURN TO APPLICATION** when ready.

Press Item Settings to begin weighing.

695.43 g

Capacity Guide: 0 to 10200.00

Diff. (wt.): -249.58 g Initial wt.: 945.01 g
Diff. (%): -26.4 % Final wt.: 695.43 g

Buttons: Item Settings, Tare, Standby, Print, Zero, Sensors, AutoCal™, Menu, More...

6. A message prompts for Item 1 to be placed on the pan again, then to press **ACCEPT** to establish final weight. After **ACCEPT** is pressed, the difference between initial and final weights appears in the Reference Fields.
7. A message prompts to start the process again, if desired.

DENSITY DETERMINATION

The balance can be used to determine an object's density. Four types of density determination can be made:

- Solids more dense than water
- Solids less dense than water
- Liquid density
- Porous material (impregnated with oil)


A density determination kit is designed to be used with the balance. Illustrations in this procedure refer to the density kit. However, you may use whatever lab apparatus that suits the requirements for density measurements. The balance software includes a built-in reference density table for water at temperatures between 50°F and 86°F. Review this entire section before attempting density measurements.

When making density measurements the material should weigh at least 10 mg on an analytical balance and 100 mg on a precision balance.

Density can be determined manually (with a key press after each step) or automatically (weights are automatically recorded). Before making density measurements, establish the Application Settings.


MEASURING THE DENSITY OF A SOLID USING WATER (DEFAULT)

1. Press the **SETUP** icon to open the density determination settings screen.
2. Confirm the following setups are selected:
 - Density type: Solid
 - Liquid type: Water
 - Porous material: Off
3. Press **DONE** to return to the density determination home screen.

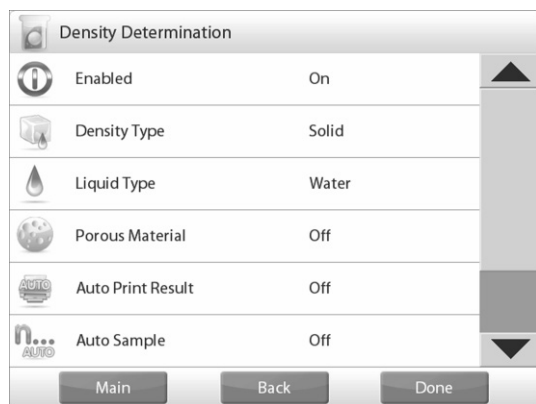EDIT SETTINGS SCREEN

- Settings available: Water temperature
 - Functions available: Return to application
1. To adjust the water temperature value, touch the **WATER TEMP** button. The balance calculates water density based on the water temperature value entered (internal lookup table). Measure the actual water temperature using a precision thermometer.

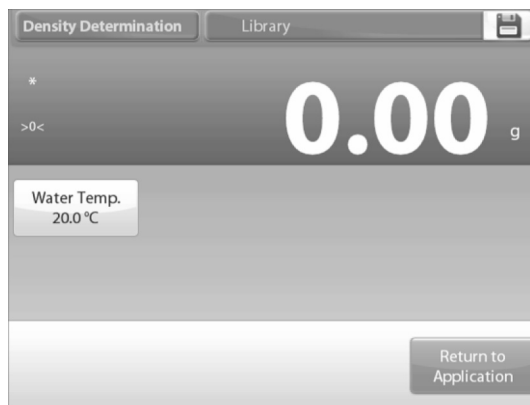2. A numeric input window appears. Enter the actual water temperature, then press **DONE**.
3. The display returns to the previous screen with the new value highlighted. To return to the density determination home screen, touch **RETURN TO APPLICATION**.
4. To begin the density determination process, press **BEGIN DENSITY DETERMINATION**.

DENSITY DETERMINATION STEPS

1. Weigh the sample in air. Press **BEGIN DENSITY CALCULATION**. Follow the screen instructions, then press **ACCEPT** to store the dry sample weight ("in air").
2. Weigh the sample submerged in the liquid. Follow the screen instructions, then press **ACCEPT** to store the sample weight (submerged in liquid).

OPERATION CONTINUED

- Once the necessary weights have been determined, the density of the sample is displayed in g/cm³ (along with the weight in air, weight in water) on the application screen.
- Press **BEGIN DENSITY DETERMINATION** to reset all the weight values and restart the process.

MEASURING THE DENSITY OF A BUOYANT MATERIAL USING WATER

- Press the **SETUP** icon to open settings screen.
- Confirm the following setups are selected:
 - Density type: Solid
 - Liquid type: Water
 - Porous material: OFF
- Press **DONE** to return to the density determination home screen.
- Follow the same procedure as solid material, except in step 2, push the sample down into the liquid until it is fully submerged.

MEASURING THE DENSITY OF A SOLID USING AN AUXILIARY LIQUID

- To enable this feature, enter the density setup menu and select the following:
 - Density type: Solid
 - Liquid type: Other
 - Porous material: OFF
- Confirm the default values displayed (liquid density, etc.) are correct.
- To edit the default values, touch **EDIT SETTINGS** and the application settings screen appears.
 - Settings available: Liquid density
 - Functions available: Return to application
- To adjust the liquid density value, touch the **LIQUID DENSITY** button.
- Measure the actual liquid density of the density setup and enter it.
- A numeric input window appears. Key in the density in g/cm³, then press √.
- The display returns to the previous screen with the new value highlighted. To return to the density determination home screen, touch **RETURN TO APPLICATION**.

- Begin the density determination process as per above.

MEASURING THE DENSITY OF A LIQUID USING A CALIBRATED SINKER (NOT SUPPLIED)

- To enable this feature, enter the density setup menu and select the following: Density type: Liquid.

 **NOTE:** When the density type is set to liquid, the liquid type and porous material selections are disabled.

- Confirm that the default value displayed (sinker volume) is correct.
- To edit the default values, touch **EDIT SETTINGS** and the Settings screen is now displayed:
 - Settings available: Sinker volume
 - Functions available: Return to application
- To adjust the sinker volume value, touch the **SINKER VOLUME** button.
- A numeric input window is displayed. Key in the desired value, then press √.
- The display returns to the previous screen with the new value highlighted.
- To return to the density determination home screen, touch **RETURN TO APPLICATION**.
- Press **BEGIN DENSITY DETERMINATION** to start the process.
- Weigh the sinker in air. Follow the screen instructions, then press **ACCEPT** to store the sinker weight ("in air").

OPERATION CONTINUED


10. Weigh the sinker submerged in the test liquid. Follow the screen instructions, then press **ACCEPT** to store the sinker weight (submerged in liquid).

11. Once the necessary weights have been determined, the density of the liquid sample is displayed in g/cm³ (along with the weight in air, weight in water) on the application screen.
12. Press **BEGIN DENSITY DETERMINATION** to reset all the weight values and restart the process.

MEASURING THE DENSITY OF POROUS MATERIAL USING OIL

1. To enable this feature, enter the density determination setup menu, and set the following:
 - Density type: Solid
 - Liquid type: Water
 - Porous material: ON
2. Confirm the default values displayed (water temperature) are correct.
3. To edit the default values, touch **EDIT SETTINGS**. The settings screen appears.
 - Settings available: **WATER TEMP, OIL DENSITY**
 - Functions available: **RETURN TO APPLICATION**

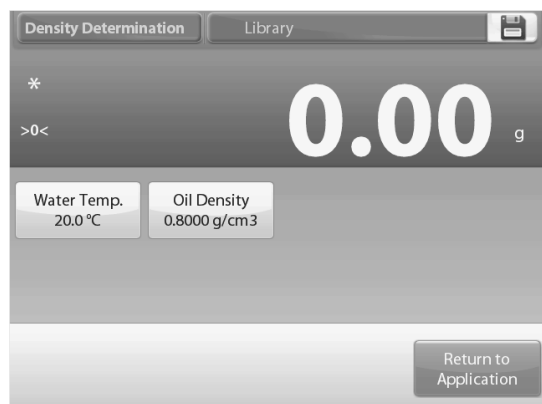4. The balance calculates water density based on the water temperature value entered (look-up table).
5. Measure the actual water temperature using a precision thermometer.
6. To adjust the water temperature or oil density values, touch the **WATER TEMP** or **OIL DENSITY** button.
7. Numeric input windows appear. Key in the desired value, then press $\checkmark$.
8. The display returns to the previous screen with the new value highlighted.
9. To return to the density determination home screen, touch **RETURN TO APPLICATION**.
10. Press **BEGIN DENSITY CALCULATION**.
11. Weigh the un-oiled sample in air. Follow the screen instructions then press **ACCEPT** to store the dry sample weight ("in air").

12. Weigh the oiled sample in air. Follow the screen instructions then press **ACCEPT** to store the sample weight (oiled).

OPERATION CONTINUED

13. Weigh the oiled sample submerged in liquid. Follow the screen instructions then press **ACCEPT** to store the oiled sample weight (submerged in liquid).

14. Once the necessary weights have been determined, the density of the sample is displayed in g/cm³ (along with the weight in air, un-oiled and oiled, and weight in water) on the application screen.
15. The value stays on the display until **BEGIN DENSITY DETERMINATION** is touched.
16. Press **BEGIN DENSITY DETERMINATION** to reset all weight values and restart the process.

PEAK HOLD

The peak hold application allows the user to capture and store the highest stable or unstable weight in a series of weight measurements.

Three different start/reset modes can be selected:

- Manual (start and stop the capture via key press)
- Semi-Automatic (start automatically/reset manually)
- Automatic (start and stop the capture automatically)

This application allows both stable and unstable weights to be captured.

PEAK HOLD – SEMI-AUTOMATIC (DEFAULT)


1. In the upper left portion of the home screen, select **PEAK HOLD**. Follow the screen instructions to begin.
2. Place item on the pan. The capture process starts automatically. The highest reading is held on the display.

3. To release the reading, press **END PEAK HOLD**.
4. Clear the pan or press **TARE** to reset the display and re-start the process.

PEAK HOLD – MANUAL

1. For manual start and reset, choose manual mode.
2. To enable this feature, enter the peak hold setup menu and change the mode to manual.

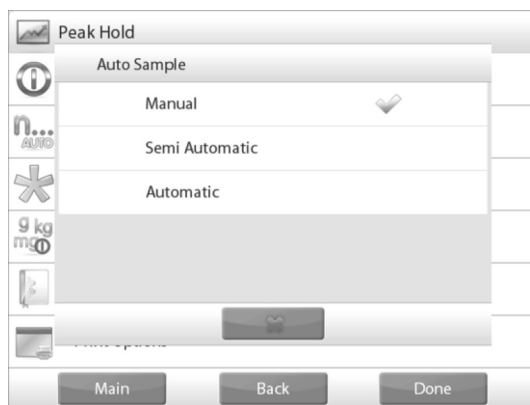OPERATION CONTINUED

THE PEAK HOLD – MANUAL HOME SCREEN


1. Follow the screen instructions: Place item on the pan.
2. Press **BEGIN PEAK HOLD** to begin the capture process.


3. The highest stable weight reading is held indefinitely.
4. Press **END PEAK HOLD** then clear the pan to reset the display and re-start the process.

PEAK HOLD – AUTOMATIC

1. For automatic start and reset, choose Automatic mode.
2. To enable this feature, enter the peak hold setup menu. Change the mode to automatic.

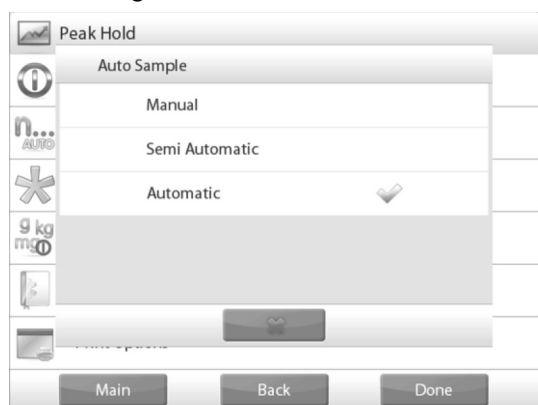3. Place item on the pan. The capture process starts automatically.

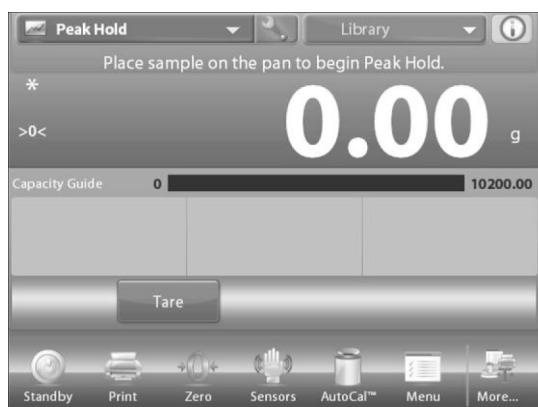4. The highest stable weight reading is held until the weight is removed.
5. The process starts automatically once the weight is removed from the pan (after 10-second delay).


INGREDIENT COSTING

This application allows the user to determine the cost of a formula or recipe based on the known cost and quantity of the components or ingredients. Up to 99 components are allowed.

1. In the upper left portion of the home screen, select **INGREDIENT COSTING**. Follow the screen instructions to begin.
2. Press **ITEM SETTINGS** to set up items in memory.

The Edit Settings Screen

- Settings available: **ADD NEW, EDIT NAME, EDIT COST, EDIT UNIT, DELETE ITEM, CLEAR ALL**
 - Functions available: **RETURN TO APPLICATION**
3. To set up the item(s) in memory, press **ADD NEW**. A new item is added (default name) with default info (cost per weight and weight unit).

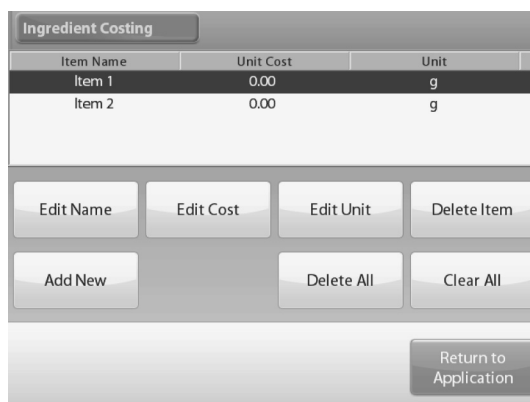4. To change the default item name, press **EDIT NAME**. The alpha-numeric entry screen appears.
5. Enter the new name, then press √. The display returns to the previous screen with the new value highlighted momentarily.
6. To change the default component cost per unit, press **EDIT COST**. The numeric entry screen appears.
7. Enter the new weight, then press √. The display returns to the previous screen with the new value highlighted momentarily.
8. To change the default component unit, press **EDIT UNIT**. The Weighting Units screen appears.
9. Continue adding items as required, using **ADD NEW**. New items can be added at any time.
10. To delete an item, touch the item to highlighting it, then press **DELETE ITEM**.

OPERATION CONTINUED

- To clear all items, press **CLEAR ALL**. A confirmation is required. When finished, press **RETURN TO APPLICATION**.
- Once setup is complete and the application screen is displayed, follow the screen instructions.
- Add the first item then press **ACCEPT** to store the value in memory. The weight of the item is recorded, then the screen prompts for the next component. Unit Cost, Total Weight and Total Cost are also displayed.

- Continue to follow the screen instructions to add the items to the pan and touch the **ACCEPT** button until each item has been weighed.
- Press the **VIEW DATA** button at any time to show the current results: Item, Item Cost and Item Weight.
- After the last item has been weighted, press **VIEW DATA** to see the results.

 **NOTE:** Cancel the current data set by pressing **CLEAR DATA**.

INGREDIENT COSTING RESULTS

- Press **PRINT** to send the data to the enabled COM port.
- Press **CLOSE** to close the results window and return to the Application home screen.


PIPETTE ADJUSTMENT

Pipette adjustment checks the accuracy and precision values of pipettes by weight analysis. The balance has a built-in density reference table for water at temperatures between 50°F and 86°F. If other liquids are used for pipette calibration, enter the liquid's density at current room temperature in g/cm³. Since all calculations are made within the balance, the atmospheric pressure is also required.

Pipette adjustment can be accomplished manually (with a key press after each step) or automatically (weights are automatically recorded when stable). Calculations made within the balance provide results of inaccuracy and imprecision. Up to 30 samples can be used.

PIPETTE ADJUSTMENT – MANUAL (DEFAULT)

- Confirm the default values displayed are correct:
 - Nominal, Inaccuracy, Imprecision, Density, Atmospheric Pressure, Water Temperature, etc.
- To edit the default values, touch the **EDIT SETTINGS** button.
 - Settings available: **NOMINAL CAPACITY, WATER TEMP, INACCURACY, ATMOSPHERIC PRESS, IMPRECISION, PIPETTE NAME, PIPETTE NUMBER**
 - Functions available: **RETURN TO APPLICATION**
- To adjust the pipette nominal capacity value, touch the **NOMINAL** button.

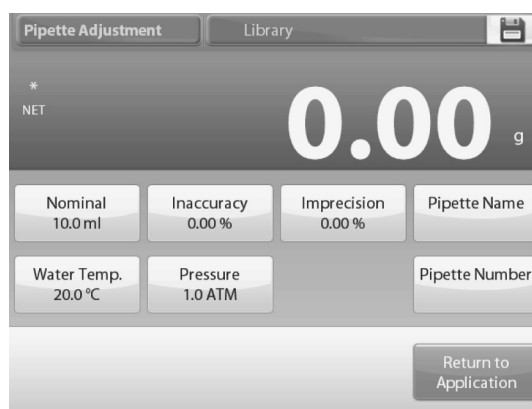
- Enter the nominal capacity in the keyboard screen that appears, and press $\sqrt{}$.
- The display returns to the previous screen with the new value highlighted.

 **NOTE:** Capacity can be measured in mL or μ L, depending on the application setup.

OPERATION CONTINUED

- To adjust the water temperature, touch the **WATER TEMP** button.
- The balance calculates water density based on the water temperature value entered. Measure the actual water temperature using a precision thermometer.
- A numeric input window is displayed.
- Key in the desired temperature, then press $\checkmark$.
- The display returns to the previous screen with the new value highlighted.
- Continue to enter the required information following the above procedure.

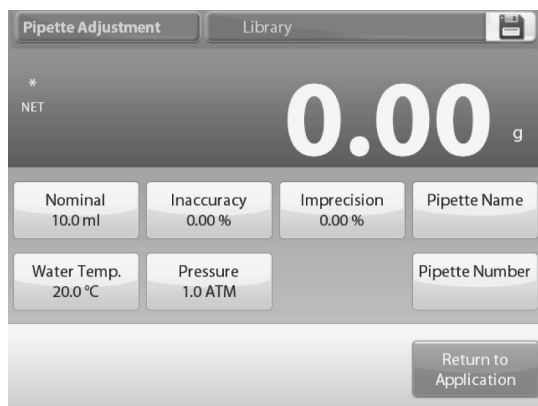- When finished, press **RETURN TO APPLICATION**. The application home screen appears, and you can begin the pipette adjustment process.

BEGIN PIPETTE ADJUSTMENT PROCESS

- Touch **BEGIN PIPETTE ADJUSTMENT**.
- Follow the screen instructions, to place container on the pan. Press **TARE**.


DISPENSE THE FIRST SAMPLE

- Follow the screen instructions, then press **ACCEPT** to store the first sample weight (liquid).
- Continue to follow the screen instructions by dispensing samples and pressing **ACCEPT** to store each sample weight (liquid).


- Once all the sample weights have been dispensed, the results screen is automatically displayed.

 **NOTE:** The default number of samples is 10 but can be changed in Setup from 2 to 30 samples.

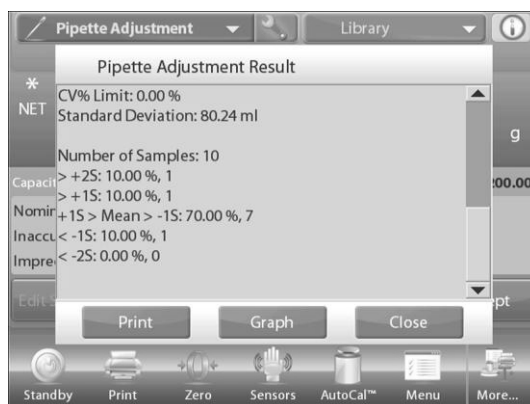- Press **GRAPH** to see results plotted. Press **RESULT** to return to the result screen. Press **CLOSE** to return to the pipette adjustment home screen.

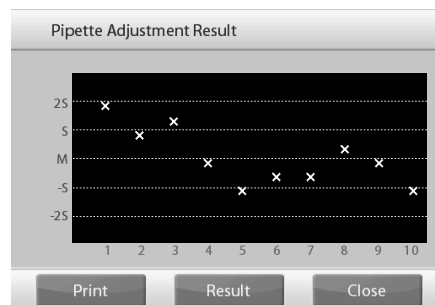- To view either the results or graph screens again from the home screen, press the **VIEW RESULT** button. To start a new pipette adjustment process, press **BEGIN PIPETTE ADJUSTMENT**.

OPERATION CONTINUED

STATISTICAL QUALITY CONTROL (SQC)

The statistical quality control (SQC) feature is useful during various types of process filling operations, to monitor and/or control the process to eliminate under- and over-filling.

As samples are weighed and stored in the balance, a trend analysis is developed and can be displayed.

Up to 10 batches with up to 80 samples in a batch are visible for quality control purposes. Each batch of samples is shown, with the maximum, minimum, standard deviation and mean values for each batch. An ongoing examination of the relative deviation of the samples along with other statistical data can be viewed.

The filling process operation can be effectively monitored by observing the results on the view trace screen.

SQC collection can be done manually (manual mode, requires key press after each step) or automatically (auto mode, the weights are automatically recorded when stable).

In the upper left portion of the home screen, select **SQC**.

EDIT SETTINGS

1. Confirm the default values (Nominal, Inaccuracy, Imprecision, Density, Atmospheric Pressure, Water Temperature, etc.) are correct.
2. To view/edit the default values, touch the **EDIT SETTINGS** button.
3. The Settings screen appears.
 - Settings Available: **NOMINAL TARGET WEIGHT**, **+TOLERANCE 1**, **-TOLERANCE 1**, **+TOLERANCE 2**, **-TOLERANCE 2**
 - Functions Available: **RETURN TO APPLICATION**

4. To adjust the Nominal Target weight value, touch the **NOMINAL** button. A numeric input window appears.
5. Key in the desired Nominal Weight, then press $\checkmark$. The display returns to the previous screen with the new value highlighted.

6. Continue to enter the required information following the above procedure:

- Positive Tolerance 1, via the **+TOLERANCE 1** button.
- Negative Tolerance 1, via the **-TOLERANCE 1** button.
- Positive Tolerance 2, via the **+TOLERANCE 2** button.
- Negative Tolerance 2, via the **-TOLERANCE 2** button.

NOTE: Tolerances can be in either weight (default) or percent.

7. When finished, press **RETURN TO APPLICATION** and the SQC home screen appears.

BEGIN THE SQC PROCESS

1. Touch the **BEGIN NEW BATCH** button.


2. First batch: the screen prompts for the batch name.
3. Enter the batch name, then press $\checkmark$.
4. The process screen is then displayed. Follow the screen instructions, then press **ACCEPT** to store the sample weight.

5. The current number of samples, along with the average, sum, minimum value, maximum value and range are displayed.
6. Continue to follow the screen instructions by weighing samples and pressing **ACCEPT** to store each sample's weight. The sample size is 10 (default). It can be changed in Setup from 2 to 80 pieces.

OPERATION CONTINUED

- At any time, press **END BATCH** to see the individual results screen, showing results for the current batch.

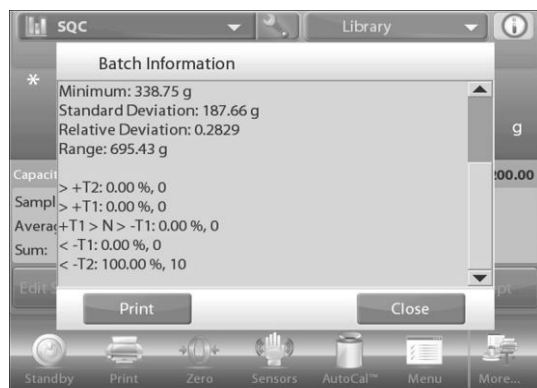
- Once all the samples have been weighed, the results screen appears automatically.
- If the COM interface is enabled, press **PRINT** to send the results either to a computer or a printer.
- Press **CLOSE** to return to the current batch process.
- A message offers the chance to save the batch, if desired.
- To view the process results of this batch (and any other saved batch), press the **VIEW BATCH HISTORY** button.

SQC				
Batch Name	Maximum	Minimum	Average	σ (stdev)
batch 1	1034.18 g	338.75 g	663.28 g	187.66 g

Buttons: Select, Clear Batches, Close

- The Batch History screen shows the basic information – maximum, minimum, average (mean), and standard deviation – for all currently saved batches.
- Press **SELECT** to see (and print) the Batch Information screen for any saved batch.

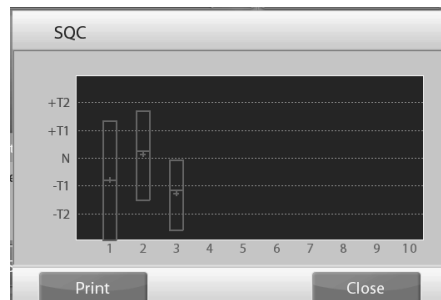- Press **CLEAR BATCHES** to reset SQC results and return to the SQC home screen.
- Press **CLOSE** to return to the SQC home screen to begin another batch or to view process results.
- To view the process results graphically (batch comparison), press the **VIEW TRACE** button.
- Press **PRINT** or **CLOSE** as shown above.

OPERATION CONTINUED

STATISTICS

Statistics are used to compare a number of samples and examine the relative deviation of the samples along with other statistical data. A minimum of three samples is required.

Statistics can operate manually (key press after each step) or automatically (weights automatically recorded when stable). To use the statistics feature, enable it in the application setup menu.

This feature is available in weighing, counting, check weighing, dynamic weighing and filling applications.

LIBRARY

When an item is processed on a regular basis, the item's data may be stored in memory for future use. This memory is referred to as the balance's library.

The following data is stored for the application used:


APPLICATION	PARAMETERS AND SETUPS	NAME	PRESET TARE	APW	REF WT. OR TARGET	SAMPLE SIZE	UNDER LIMIT	OVER LIMIT	LEVEL	MODE
Weighing	x	x	x							
Counting	x	x	x	x		x	x	x		x
Percent Weighing	x	x	x		x					
Check Weighing	x	x	x				x	x		x
Dynamic Weighing		x	x						x	x
Filling	x	x	x		x					
Totalization	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Formulation	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Differential Weighing	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Density Determination	x	x	x							
Peak Hold	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Ingredient Costing	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Pipette Adjustment	x	x	x							
SQC	x	x	x							

CREATING A LIBRARY RECORD

1. To create a library record, the settings screen is used. From the home screen, touch the **EDIT SETTINGS** button. The settings screen is now displayed.
2. Touch the **EDIT SETTINGS** button.
3. The Settings screen appears. Confirm the application settings; then touch **SAVE** icon.
4. An alpha-numeric entry screen appears to enter the library record name.
5. Load a library record from the home screen.
6. Enter the record name, then press √.

NOTE: Library record names can be eight characters or less.

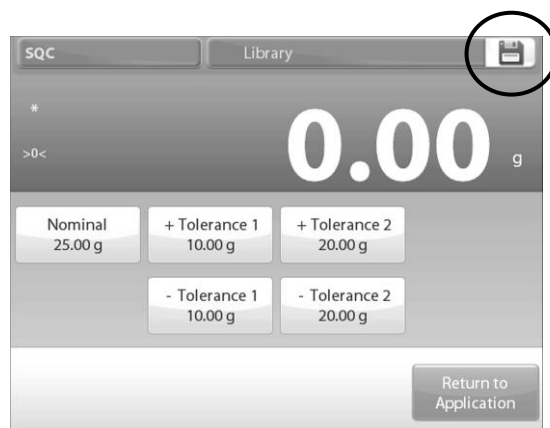7. A record is saved message appears. Press **OK** then press **RETURN TO APPLICATION**.

OPERATION CONTINUED

RETRIEVING A LIBRARY RECORD

Load a library record from the home screen.

1. Press the **LIBRARY** bar in the upper right area of the screen. It shows the name of the most recently used library file. The library manager pop-up is now displayed.

2. Touch the library record to highlight it, then touch **RECALL** to load the record's settings into the application.
3. Then press **CLOSE** to return to the home screen with the library name displayed and the stored parameters active.

DELETING A STORED LIBRARY RECORD

1. To delete a stored record, follow Retrieving a Library Record above, but press **DELETE** instead of **RECALL** in step 2.

WEIGH BELOW

The balance is equipped with a weigh below hook for weighing below the balance.


NOTE: Before turning the balance over, remove the pan, pan support, wind ring/shield and draft shield elements (if present) to prevent damage.

CAUTION! Do not place the balance on the pan support cone or load cell pins.

To use this feature, remove power from the balance, then remove the protective cover for the weigh below opening (two screws). The protective cover is reversible for easy storage.

The balance can be supported using lab jacks or any other convenient method. Ensure the balance is level and secure. Power on the balance, then use a string or wire to attach items to be weighed.

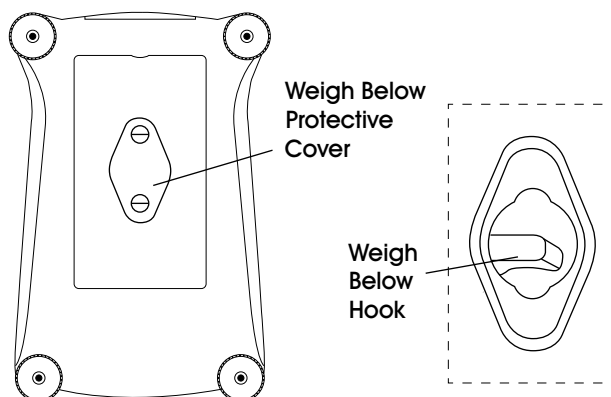MAINTENANCE

CALIBRATION

Periodically verify calibration by placing an accurate weight on the balance and viewing the result. If calibration is required, perform a balance internal calibration.

CLEANING

CAUTION! Disconnect the balance from the power supply before cleaning. Ensure no liquid enters the interior of the terminal or base.

1. Clean the balance at regular intervals.
2. Housing surfaces may be cleaned with a lint-free cloth slightly dampened with water or a mild cleaning agent.
3. Glass surfaces may be cleaned with a commercial glass cleaner.
4. Do not use solvents, harsh chemicals, ammonia or abrasive cleaning agents.


TROUBLESHOOTING

OPERATING ISSUES	POSSIBLE CAUSES	RECOMMENDATIONS
Balance will not turn on.	No power to balance.	Verify connection and voltage.
Poor accuracy.	Improper calibration. Unstable environment.	Perform calibration. Move balance to suitable location.
Cannot calibrate.	Calibration menu is locked. Unstable environment. Incorrect calibration masses.	Turn calibration menu lock off. Move balance to suitable location. Use correct calibration masses.
Cannot change menu settings.	Sub-menu is locked.	Unlock sub-menu.
Low reference weight.	Reference weight too small. The weight on the pan is too small to define a valid reference weight.	Increase sample size.
Invalid piece weight.	Average piece weight is too small.	Increase average piece weight.
Operation time out.	Weight reading is not stable.	Move balance to suitable location.
Screen displays "-----".	Busy (tare, zero, printing).	Wait until completion.
Scale will not weigh or is displaying inaccurate weights.	Scale has fallen out of calibration.	Recalibrate scale. Refer to page 5.

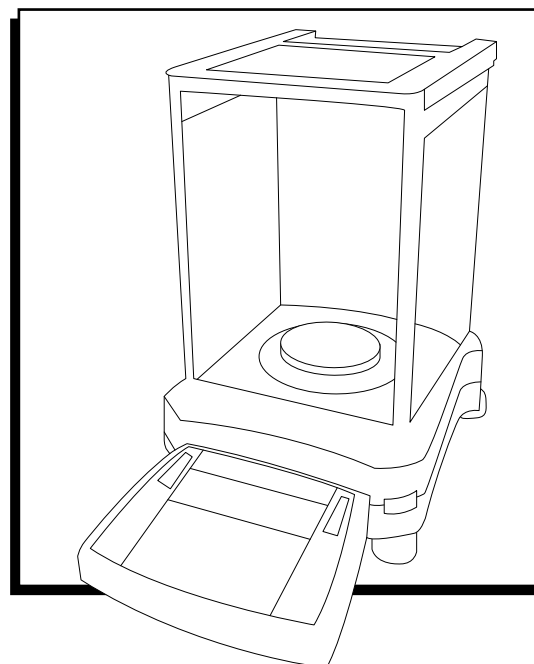If the troubleshooting section does not resolve your problem, contact Uline Customer Service at 1-800-295-5510.

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