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
1-800-295-5510
uline.com

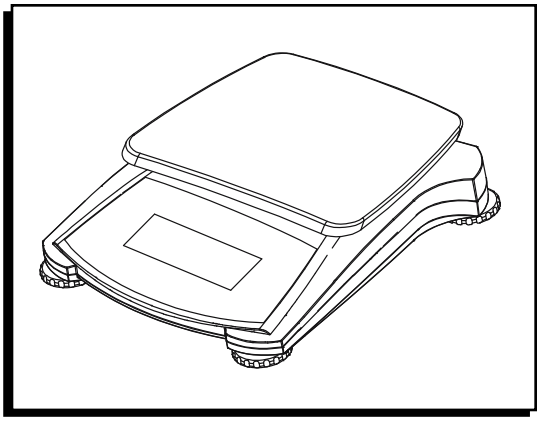TROUBLESHOOTING


OPERATING ISSUES	RECOMMENDATIONS
Cannot turn on	No power to scale.
Poor accuracy	Improper calibration. Unstable environment.
Not displaying two decimal places	Security switch in locked position. Move to unlocked position.
Cannot calibrate	Security switch in locked position.
Not displaying two decimal places.	Security switch in locked position. Move to unlocked position.
Scale will not weigh or is displaying inaccurate weights.	Transportation lock is in locked position. Move to unlocked position. Scale may have fallen out of calibration. Recalibrate scale. Refer to page 6 of instructions.

ERROR CODES	RECOMMENDATIONS
Err 8.1	Weight exceeds power on zero range.
Err 8.2	Weight below power on zero range.
Err 8.3	Overload (weight exceeds rated capacity).
Err 8.4	Underload (pan removed).
Err 8.5	Tare out of range.
Err 8.6	Displayed value exceeds 999999 (possible in Totalization mode).
Err 8.7	Unknown internal weight position.
rEF.Err	Parts counting or percentage error – sample weight <1d. Scale shows error, then exits parts counting or goes to [Clr.APU].
Lo.rEF	Percent Reference weight or APW is too low for accurate results.
CAL E	Failure to calibrate.
USb.Err	Cannot find menu or app file in U-disk.

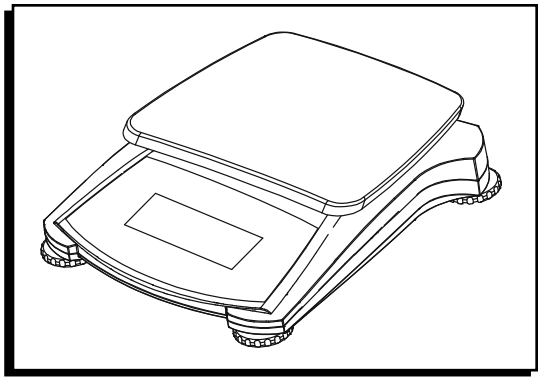If the troubleshooting section does not resolve or describe problem, contact Uline Customer Service at 1-800-295-5510.

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